

UNITED SKIN CARE

QUANTUS

LEO

- 1200 Watt 4 Wave Diode Laser
- 1000 W Qswitch Ndyag
- 1 Year Warranty
- Ce Approved
- Lhr Results In 6-7 Sessions
- Pan India Service
- 5 Wavelength Pico
- Medical Ce Approved



Triple Wave
Diode



NDYAG
Laser



MACHINE DISPLAY



Therapy Principle Diode

Diode laser hair removal System, adopts the light selective absorption theory, heating the hair shaft and hair follicle to destroy the hair follicle and oxygen organization around hair follicle, thereby to get the permanent hair removal purpose.

The patent Pulse Width to Energy 1:1 working mode, makes the energy penetrating more effective, less thermal destroy to the surrounding tissue, meanwhile, the stable anti-ice cooling technology, protects skin from burning finally getting a wonderful result.



Therapy Principle Pico

THERAPY OF EPIDERMIS AND DERMIS PIGMENTATION

Utilizing the explosive effect of the Nd: YAG laser, the laser permeate the epidermis into the dermis which includes amount of pigment mass. Since the laser pulses in nanosecond but with super high energy, the shot pigment mass swell quickly and breaks into small pieces, which will be eliminated through metabolic system.

CARBON PEELING TREATMENT

The 1320 nm wavelength laser directly acts on the carbon cream, blasting the toner, and the toner absorbs pore oil and dirt. At the samea time, the laser can stimulate the dermis to produce a large amount of collagen, thereby achieving the function of shrinking pores and whitening.



Treatment Range

- Skin Rejuvenation, Skin Whitening and Resurfacing.
- Freckles, Black spot, Age spots sun-burn spots, Coffee spots.
- Pigment of epidermis and derma layer.
- Body tattoo, eyebrow tattoo and lips line.
- Nevus of Ota and birthmark.

Pico Laser Filters



1064 nm Head

Black, blue, deeper penetration
Melasma



532 nm Head

532nm: Red, orange, brown,
stronger exfoliation.



Carbon facial Head

Carbon Peeling Head



1064 nm lite head

Eyebrow Tattoo Removal



532 nm lite Head

Lip Tattoo Removal

System Advantages



Precise

The laser generator is located at the handle position, which is not prone to polarization.



Easy To Use

Simple Software easy to operate.



High Performance

Short pulse width(4000ps) High power(2000mj) Long lifespan



Handle Design

Sleek and light weight handle



Good Cooling

High power built in water cooling system make the machine reliable for continually use.



Color Display

Operation and parameter adjustment with 10.4 inch color TFT touch screen.

- **LASER TYPE** : 4 Diode + ND: yag Laser passive
- **POWER** : 2500W
- **OUTPUT SYSTEM LASER** : 1200 w Handle tool output
- **WAVELENGTH** : 4 Wave Diode and 1064nm,532nm,and 1320nm
- **PULSE WIDTH** : 5ns
- **ENERGY** : 100-2000mj
- **FREQUENCY** : 1-10hz
- **COOLANT** : Distilled water
- **SCREEN** : 15.6 inch touch color screen
- **COOLING SYSTEM** : Air+water cooling
- **VOLTAGE** : 110V±10%, 60Hz /220V±10%, 50Hz
- **MACHINE WEIGHT** : 50kg



**Female | Age
26 Before**

**After
5 treatments
2019.08.10-
2020.01.20**

**Male | Age
39 Before**

**After
4 treatments
2019.08.08-
2020.01.01**





Before




After



Before



After

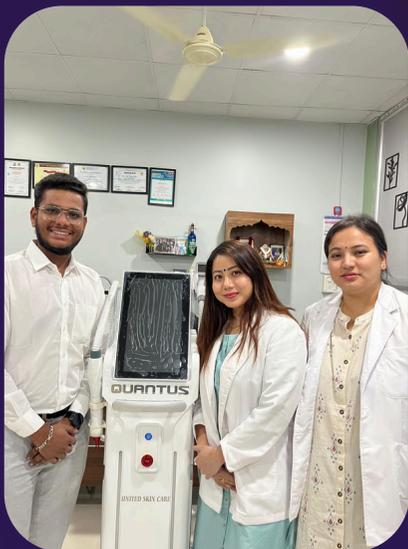

UNITED SKIN CARE

PAN INDIA INSTALLATIONS



UNITED SKIN CARE



Onsite Training

Professional training of clinical operations and technical maintenance

1 Year Warranty

Onsite 1 year Warranty And pan India Service



UNITED SKIN CARE

TO ORDER CALL OR WHATSAPP ON
8160903565/9284416452